

IGRT kV-imaging dose MC calculations validated in anthropomorphic phantoms using OSL

L. Berger, Guillaume Boissonnat, Helena Chesneau, V. Passal, J. Desrousseaux, S. Gempp Gempp Raucoules, G. Delpon, O. Henry, J. Bellec, Thomas Dautremer, et al.

▶ To cite this version:

L. Berger, Guillaume Boissonnat, Helena Chesneau, V. Passal, J. Desrousseaux, et al.. IGRT kV-imaging dose MC calculations validated in anthropomorphic phantoms using OSL. ESTRO 37 - European SocieTy for Radiotherapy & Oncology, Apr 2018, Barcelona, Spain. Radiotherapy and Oncology, 127 (Supplement 1), S1004-S1005; poster EP-1860, 2018, ESTRO 37, April 20-24, 2018, Barcelona, Spain. 10.1016/S0167-8140(18)32169-8. hal-01878965

HAL Id: hal-01878965 https://univ-rennes.hal.science/hal-01878965

Submitted on 21 Aug 2019

HAL is a multi-disciplinary open access archive for the deposit and dissemination of scientific research documents, whether they are published or not. The documents may come from teaching and research institutions in France or abroad, or from public or private research centers.

L'archive ouverte pluridisciplinaire **HAL**, est destinée au dépôt et à la diffusion de documents scientifiques de niveau recherche, publiés ou non, émanant des établissements d'enseignement et de recherche français ou étrangers, des laboratoires publics ou privés.





IGRT kV-CBCT dose MC calculations validated in anthropomorphic phantoms using OSL





cea

L. BERGER¹, G. Boissonnat², H. Chesneau³, V. Passal⁴, J. Desrousseaux⁵, S. Gempp Raucoules⁵, G. Delpon⁴, O. Henry³, J. Bellec³, T. Dautremer², E. Barat², J. Garcia-Hernandez², L. Padovani⁵, C. Lafond³, D. Lazaro²


¹Centre Jean Perrin, Radiotherapy, CLERMONT FERRAND, France; ²CEA, LIST, F-91191 GIF SUR YVETTE, France; ³Centre Eugène Marquis, Inserm U1099- LTSI, Rennes, France; ⁴Centre René-Gauducheau- Institut de Cancérologie de l'Ouest, Radiotherapy, Saint-Herblain, France ; ⁵Assistance Publique - Hôpitaux de Marseille, Radiotherapy, Marseille, France.


While in-room Magnetic Resonance Imaging starts becoming part of radiotherapy (RT) treatments, the use of X-ray imaging equipment in Image-Guided RT (IGRT) is still growing and with it, the need to evaluate the additional dose-to-organs it delivers. This study aims to verify the accuracy of Monte Carlo (MC) calculation of the patient's dose-to-organs delivered by four commercially available kV imaging systems: XVI CBCT (Elekta), OBI CBCT (Varian), ExacTrac 2D-kV system (Brainlab) and 2D-kV CyberKnife imaging system (Accuray). Simulations were validated against OSL measurements in the pediatric anthropomorphic phantom Grant (CIRS, ATOM) performed in three different clinical sites.


kV X-ray source description

All four systems kV X-ray source were described following the same principles:

- **Stage 1:** Geometry implementation and MC calculation (Penelope)
- ⇒ Electron incoming energy and direction
- ⇒ Tungsten anode angle
- ⇒ Initial filtration, additional filtration (bowtie filters)
- ⇒ With no physical collimation
- ⇒ The **photon energy fluence distribution** is finally stored in a phase space file further away from the origin
- **Stage 2: Photon energy fluence distribution** simplification
- ⇒ Photons geometrically outside of the virtual selected collimation are removed from the distribution.
- \Rightarrow The photon energy fluence distribution $\phi(x, y, E)$ is simplified as the product of two 1-D spatial functions $p_x(x)$ and $p_y(y)$ and a 2-D spatial-energetical function $p_{E|x \text{ or } y}(x \text{ or } y, E)$
- ⇒ The three functions are smoothed using taking into account the MC noise
- **Stage 3:** Generating new photons from the virtual source model
- ⇒ The cumulative functions of the three smoothed distributions are used to sample the position and energy of new photons. ⇒ Obtained photons are translated and rotated according to
- every X-ray tube incidence used by the simulated imaging protocol.

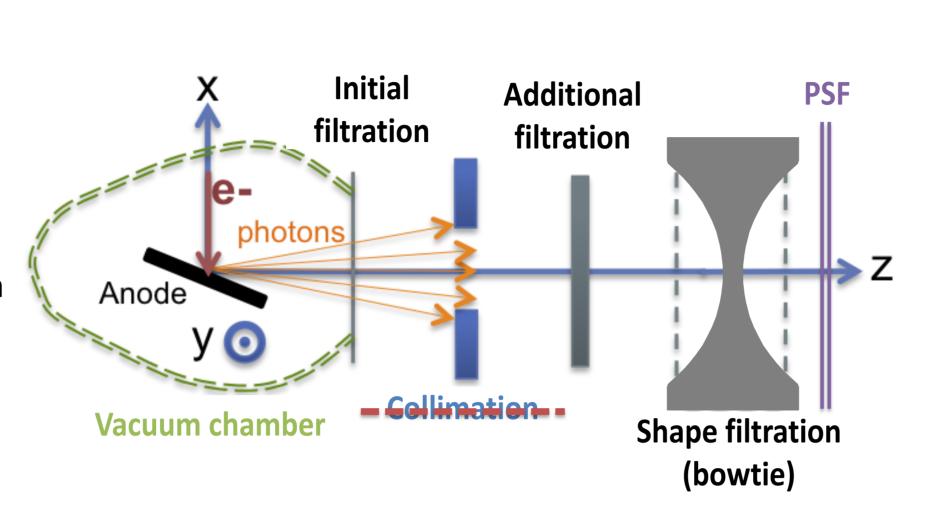

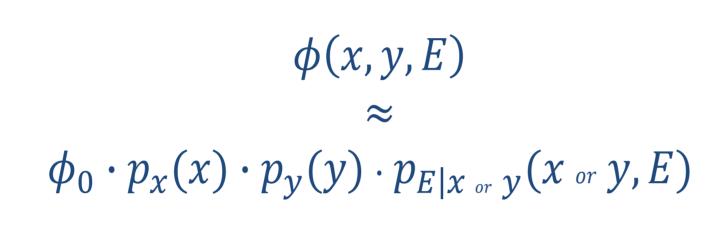

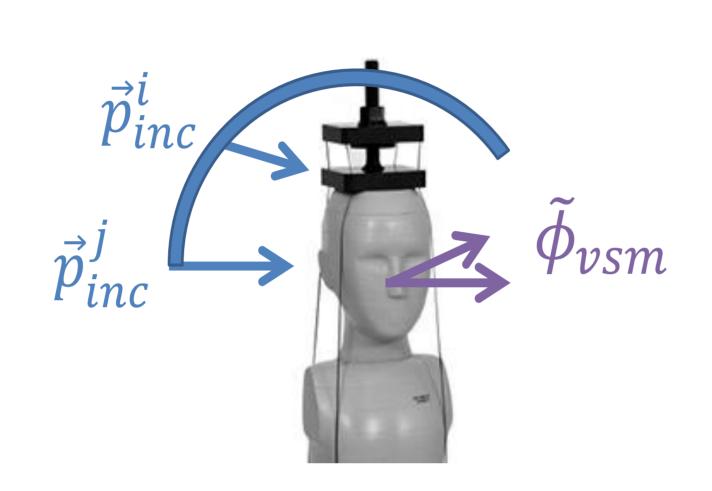



MC calculations comparison with OSL measurements


This study includes OSL measurements performed in radiotherapy centers with pediatric protocols* for the four IGRT. Three exam localizations were used: head (15 OSL), thorax (17 OSL) and pelvis (13 OSL). Monte Carlo dose maps were calculated in 2 hours on 40 CPU using an in-house parallelized version of Penelope.


Taking in account measurements, calculations and positioning uncertainties, agreement between OSL measurements and MC calculations a 10%-3mm Gamma Index criteria was used.

Over all the exams in the study, the Gamma index has a 87% success rate (94% if removing points localized outside of the geometrical collimation).

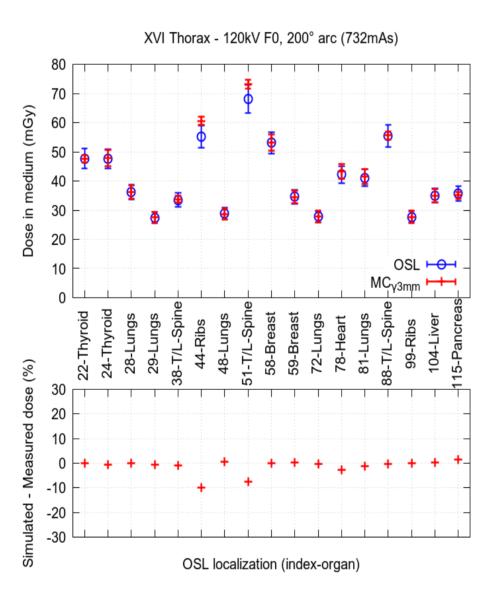

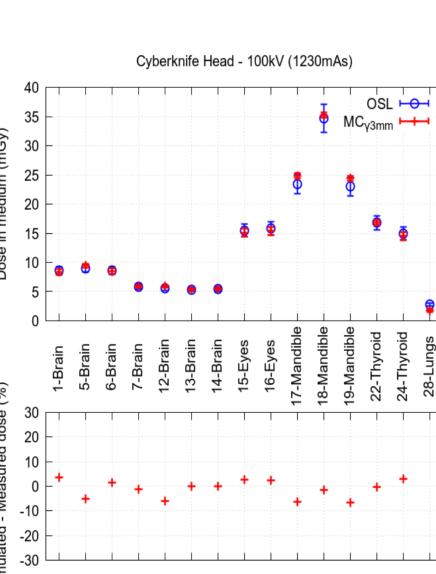

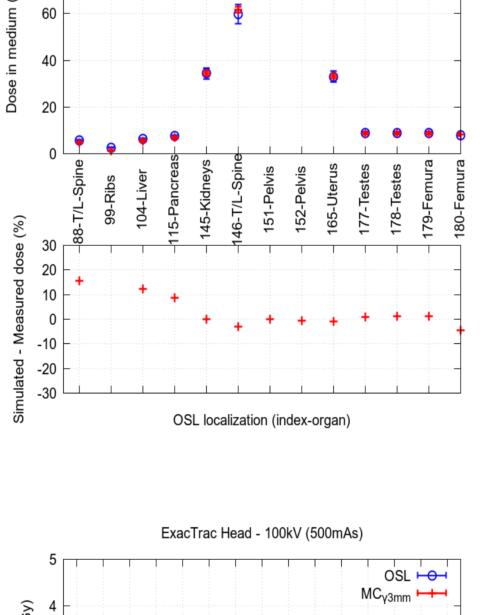

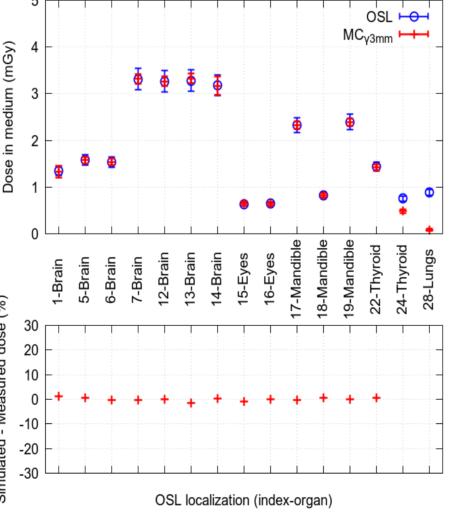

DOI: 10.3252/pso.eu.ESTRO37.2018


*beam intensity (mAs) were increased compared to standard protocols in order to decrease OSL dose measurement uncertainties

Voxelized anthropomorphic phantom construction


The Pediatric anthropomorphic phantom: Grant (CIRS, ATOM) was scanned. The obtained 3D HU map was then converted into a 3D material map using the chemical composition found in the Grant documentation.





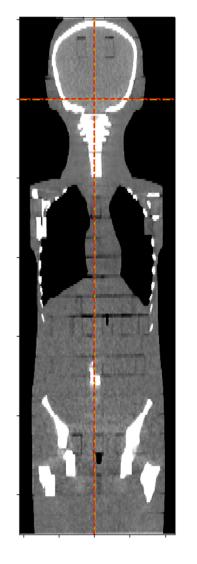

Material assignation

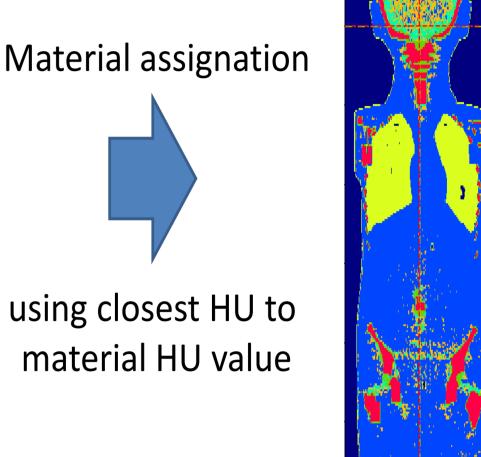


The in-house Penelope version (CEA-LIST) allowed us to simulate the dose deposition using the virtual source model in this Grant voxelized model.

OSL measurement


In this study the OSL nanodots detectors from Laundauer were used.

- **Step 1: OSL air-Kerma cross calibration**
- \Rightarrow For every beam quality Q, OSL nanodots were calibrated with an air ionization chamber (with a 3% precision).
- **Step 2: From air-Kerma to dose in medium**
- ⇒ For comparison purposes, air-Kerma OSL measurements were converted into dose-in-medium using the conversion factor
- \Rightarrow As it is highly difficult to know the beam quality Q' at the measurement point, we assumed Q' = Q. The induced error $\varepsilon_{o'}^Q$ was estimated, through MC calculation, as being inferior to 3% (1% in average).

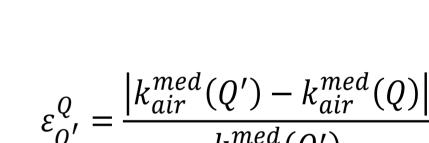

 $k_{air}^{med}(Q') =$

Dose level comparisons between systems

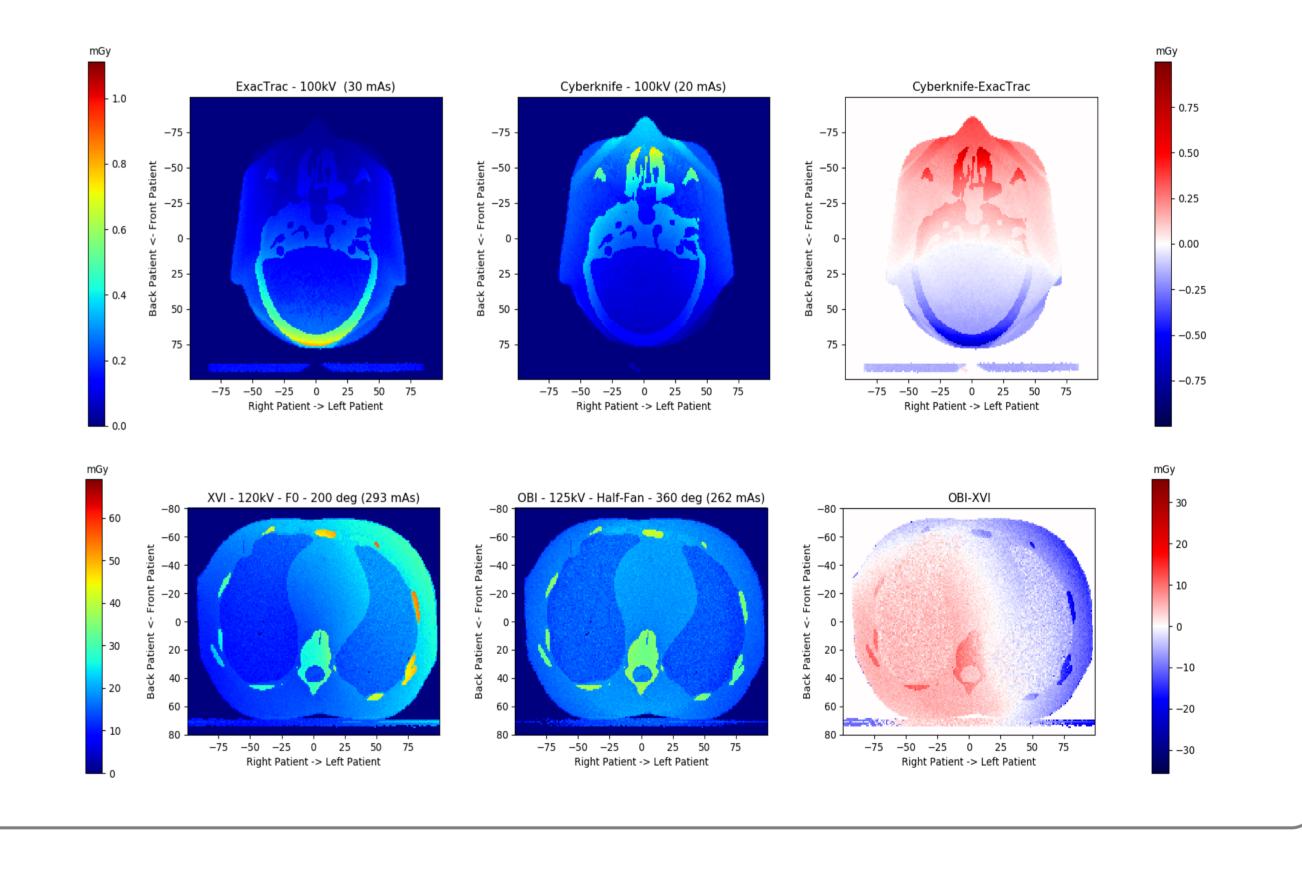

Perspectives

As part of the AID-IGRT project, an in-house software was developed which allows to:

- Transform a CT image into a voxelized material map of the patient
- Select the imaging protocol (system, voltage, collimation, filtration, irradiation arc)
- Calculate the dose-in-tissue using the Penelope MC platform

In addition, a model of the Tomotherapy MV-CT was added to the software, thanks to developments done at CRG (Nantes).

This tool is currently being tested in clinical centers in order to evaluate the dose level currently added by the use of IGRT during the course of a radiotherapy treatment.

Summary

MC dose Calculation

- The photon distributions of four kV-IGRT systems were obtained using MC and converted into Virtual Source Model
- A voxelized material map of the anthropomorphic phantom was obtained by converting HU values from a CT scan into materials
- Photons from VSM are propagated into the voxelized geometry and the corresponding dose-in-medium are calculated using Penelope

Experimental validation

- In phantom, doses were measured using OSL detectors calibrated in air Kerma for the four kV-systems
- Measured doses were converted into dose-in-medium
- A 87% Gamma Index (10%, 3mm) success rate was obtained between measured and MC-calculated doses.

Overall, the method presented here seems suitable to reproduce doses with an accuracy better than 20% for various IGRT protocols and could be used to track the impact of the current IGRT practices on RT additional doses.

Reterences

V. Passal et al. "Monte Carlo calculation of absorbed doses due to imaging sessions delivered to patients

during Tomotherapy Image-Guided RadioTherapy courses" MCMA2017, Napoli 2017. G. Boissonnat et al. "IGRT kV-CBCT dose calculations using Virtual Source

Models and validated in phantoms using OSL", ESTRO annual meeting 36, Vienna 2017.

H. Chesneau et al. "Comprehensive validation of a Monte Carlo kV-CBCT model using OSL and spectral measurements" ESTRO annual meeting 35, Torino 2016.

H. Chesneau et al. "Monte Carlo simulation for imaging dose estimation: application to the Elekta XVI kV-CBCT" ESTRO annual meeting 34, Barcelona 2015.

ACKNOWLEDGEMENTS

This work is part of the French National Research Agency (ANR) project AID-IGRT involving the research institute CEA (Saclay), the Dosisoft company (Cachan) and four French cancer treatment centers: CJP (Clermont-Ferrand), CEM (Rennes), AP-HM (Marseille) and CRG (Nantes).

